


PCS[®]2

Your choice for Plasma Collection

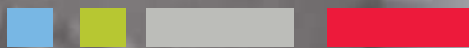

HAEMONETICS[®]
THE Blood Management Company

PCS[®]2 Plasma Collection System

Easy to use

- Self-loading pumps
- Auto-priming
- Comprehensive error and diagnostic messages
- Advanced optical sensors
- Easy-to-read display
- Minimal alarms

Multiple Protocols



- Platelet Poor Plasma
- Fresh Frozen Plasma
- All protocols with optional saline compensation

Multiple Safety Features

- Redundant air detectors
- Donor line pressure monitor
- Low AC ratios
- Fluid sensor
- Enclosed centrifuge well

Compact

- Small
- Easy to move
- Lightweight
- Ideal for centers with restricted space and for mobile collections

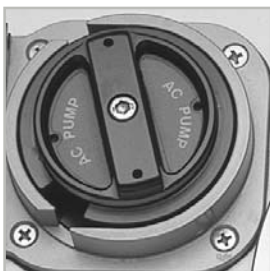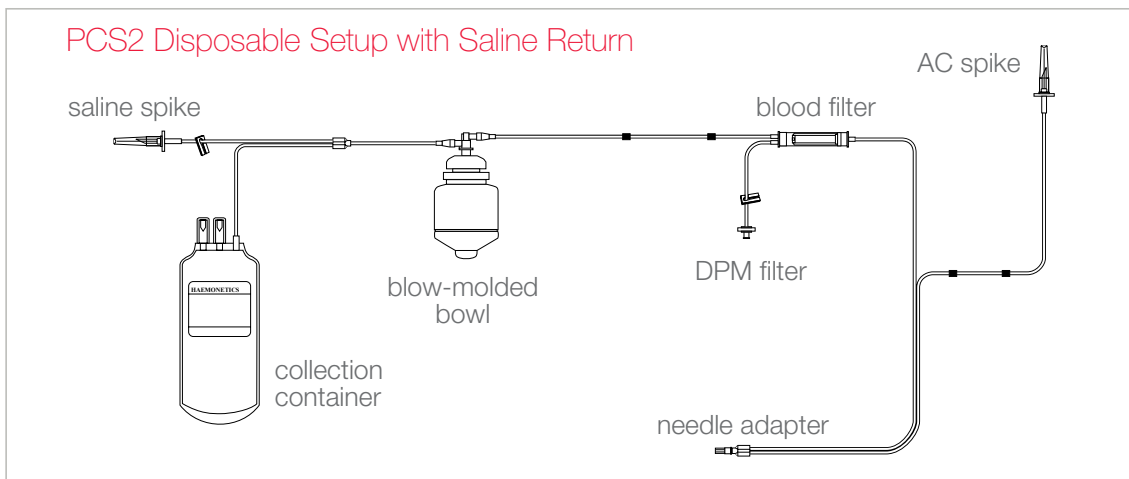

PCS[®]2 is now available with EXPRESS[™]: Intelligence and speed without compromise

FDA cleared and CE marked

What is EXPRESS?

- EXPRESS is a new software development for the PCS[®]2
- EXPRESS enhancements decrease collection time by an average of 20%, or approximately 10 minutes, per procedure¹
- EXPRESS software improves donor comfort
- EXPRESS software delivers the same or improved quality plasma product



What are the EXPRESS enhancements?

- Intelligent algorithm
 - calculates **actual** donor vein pressure, thereby optimizing the flow rate throughout the procedure
 - increases draw and return
 - improve donor experience

Why should I install EXPRESS?

- Improve donor experience by decreasing donor on-device time while using 16- and 17-gauge needles.
 - Using a smaller needle doesn't compromise the performance enhancements
 - Improved donor experience translates into improved donor retention
- Improve center production throughput by increasing number of collections during normal operating hours
 - Alternatively, may reduce center operating hours without compromising current collection volume
- Improved device management
 - Less need for center expansion
 - Reduce operator workload
 - Fewer devices to maintain or service

How do I change to EXPRESS?

Simple installation process done in center by a Haemonetics Field technician.

Will I need to retrain my staff?

No additional training is required. EXPRESS[™] enhancements are transparent upon new software installation and don't require operator intervention.

¹ Toby Simon, MD, ZLB Plasma, ASFA 2008. May not be representative of all donations.

Specification and Ordering Information for PCS2

Specifications		
PCS2 alone (closed)	H x W x D	44 cm x 55 cm x 33 cm
PCS2 with cart (open)	H x W x D	140 cm x 57 cm x 67 cm
	pole height	67 cm from top deck
Weight of PCS2		25 kg
Weight of cart		13,2 kg
Electrical voltage		230 VAC +/- 10%
Operating frequency		47-63 Hz
Current		F2,5 A @ 250 V
Power cord length		4,88 m

Description	List No.	Units / Case
PCS2 device	6002-CP-220-EW	1
Donor Blood Line	620	100
Blow-molded Bowl	625B	25
Collection Bag	690	102
Collection Bag with Vacuum Sampleadapter	SC690	70
Collection Bag with Saline Adapter	692	72

For further information regarding additional disposables, please contact your local sales representative.

The EXPRESS new Software

When ordering EXPRESS, you will receive from Haemonetics:
Per device
(1) Express program chip
(1) Express device decal
Per site
Detailed upgrade instructions
Operator's Manual
Express Launch package containing:
<ul style="list-style-type: none"> ▪ Wall poster ▪ Express documentation
Contact your local account manager for upgrade details and pricing information

Austria
0800 29 2777

Belgium and Luxembourg
(FR): 0800 754 80
(NL): 0800 754 82

Czech Republic
800 143 243

Denmark
8088 7112

France
0800 90 11 58

Germany
0800 180 8890

Italy
800 870 200

Norway
800 18 453

Sweden
020 797 150

Switzerland
0800 898 898

The Netherlands
0800 0222 707

United Kingdom
0808 2344817 or
0808 101 1375